

GENOTYPING BY PCR PROTOCOL

MUTANT MOUSE RESOURCE & RESEARCH CENTER: UC DAVIS

mmrrc@ucdavis.edu
530-754-MMRRC

NAME OF PCR: STOCK Tg(Sox3-EGFP)MT59Gsat/Mmucd MMRRC # 036040-UCD

Protocol:

Reagent/ Constituent	Volume (μL)
Water	11.275
10x Buffer	2.5
MgCl ₂ (stock concentration is 25mM)	1.7
Betaine (stock concentration is 5M)	6.5
dNTPs (stock concentration is 10mM)	0.5
DMSO	0.325
Primer 1 (stock concentration is 20μM)	0.5
Primer 2 (stock concentration is 20μM)	0.5
Taq Polymerase (5Units/μL)	0.2
DNA extracted with ☐ NaOH ☒ Proteinase K ☐ Other:	1.0
TOTAL VOLUME OF REACTION:	25μL

Comments on protocol:

- Use Touch-Down cycling protocol-first 10 cycles anneal at 65° C decreasing in temperature by 1.0° C; next 30 cycles anneal at 55° C.
- Betaine/DMSO is standardized due to high GC content in promoter regions. Also, may adjust MgCl₂ to increase reaction or decrease non specific amplifications.

Strategy:

Steps	Temp (°C)	Time (m:ss)	# of Cycles
1. Initiation/Melting HOT START? ☐	94	5:00	1
2. Denaturation	94	0:15	} 40x
3. Annealing	65 to 55 (↓1°C/cycle)	0:30	
4. Elongation		0:40	
5. Amplification	72	5:00	1
6. Finish	4	Hold	n/a

Primers:


Name	Nucleotide Sequence (5' - 3')
1: Sox3 (36040) F	GAGTCCCCCGGCCTTTTCAC
2: GFP R3	GGTCGGGGTAGCGGCTGAA

Electrophoresis Protocol:

Agarose: 1.5% V: 90

Estimated Running Time: 90 min.

Primer Combination	Band	Genotype
1 and 2	350 bp	transgenic

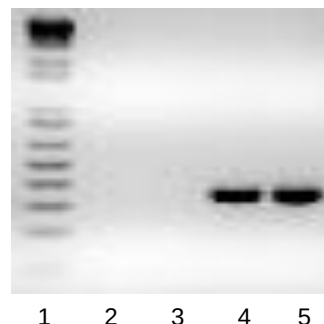

Lanes
1. 1 kb+ ladder
(Invitrogen, Cat. #10787-026)
2. Non-template control
3. Wild-type & eGFP
4. Sox3 tg/+