

KOMP PCR Design

Mouse PCR Protocol (version 1JE)

Design ID: 232256

Project ID: CSD36674

Selection Cassette: L1L2_Bact_P



MMRRC Stock #: 046395-UCD

C57BL/6N-Atm1Brd


A230050P20Riktm1.1(KOMP)Mbp/Mmucd

Suggested DNA Prep: DNeasy®Tissue Kit

Reagent	1X (μL)
water (biological grade)	10.725
betain 5M (Sigma)	6.5
DMSO (Sigma)	0.325
10X buffer w/o MgCl ₂ (AB)	2.5
25 mM MgCl ₂ (AB)	1.75
10 mM dNTPs (Invitrogen)	0.5
primers (20 μM each)	0.5
Taq 5U/μL (AmpliTaq, AB)	0.2
total cocktail	23
template	2
reaction volume	25

Cycling Parameters	
Temperature °C	Time
94	5 min
94	15 sec
65	30 sec
72	40 sec
10X (decrease 1°C/cycle)	
94	15 sec
55	30 sec
72	40 sec
30X	
72	5 min
4	finished

Primer Strategy

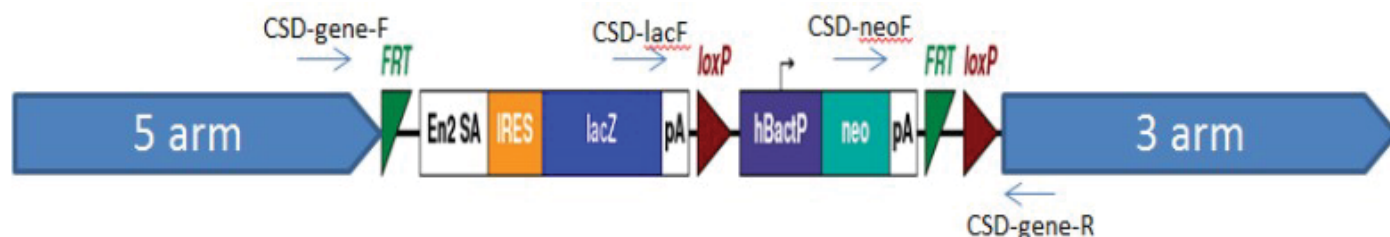

Cassette Primers

CSD-lacF: GCTACCATTACCAAGTTGGTCTGGTGTGTC
CSD-neoF: GGGATCTCATGCTGGAGTTCTTCG

Gene Specific Primers

CSD-R: GAACCAAAGGACAGAGAGATCAGAGG
CSD-F: TAGGGTCAGAGGTGAGAAGGAAAGG
CSD-wtR: GTGTGGTCGCTGCCAAACTTCC

Genotype	Forward Primer	Reverse Primer	Amplicon size (bp)
PreCre	CSD-neoF	CSD-R	632
PostCre	CSD-lacF	CSD-R	583
Wildtype	CSD-F	CSD-wtR	392
PostFlp	CSD-F	CSD-R	630

Please note, these primers are auto-designed and may not have been verified by the repository, and as such may require optimization or redesign by your facility.

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the PreCre primers may be used to identify the mutant mice. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line. The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).