

KOMP PCR Design

Mouse PCR Protocol (version 1JE)

Design ID: 88443

Project ID: CSD24480

Selection Cassette: L1L2_Bact_P



MMRRC Stock #: 046799-UCD

C57BL/6N-Bhlhe40tm1a(KOMP)Wtsi/JMmucd

Suggested DNA Prep: DNeasy®Tissue Kit


| Reagent | 1X (µL) |
|---------------------------------------|---------|
| water (biological grade) | 10.725 |
| betain 5M (Sigma) | 6.5 |
| DMSO (Sigma) | 0.325 |
| 10X buffer w/o MgCl ₂ (AB) | 2.5 |
| 25 mM MgCl ₂ (AB) | 1.75 |
| 10 mM dNTPs (Invitrogen) | 0.5 |
| primers (20 µM each) | 0.5 |
| Taq 5U/µL (AmpliTaq, AB) | 0.2 |
| total cocktail | 23 |
| template | 2 |
| reaction volume | 25 |

| Cycling Parameters | |
|--------------------|----------|
| Temperature °C | Time |
| 94 | 5 min |
| 94 | 15 sec |
| 65 | 30 sec |
| 72 | 40 sec |
| 94 | 15 sec |
| 55 | 30 sec |
| 72 | 40 sec |
| 72 | 5 min |
| 4 | finished |

10X (decrease 1°C/cycle)

30X

Primer Strategy

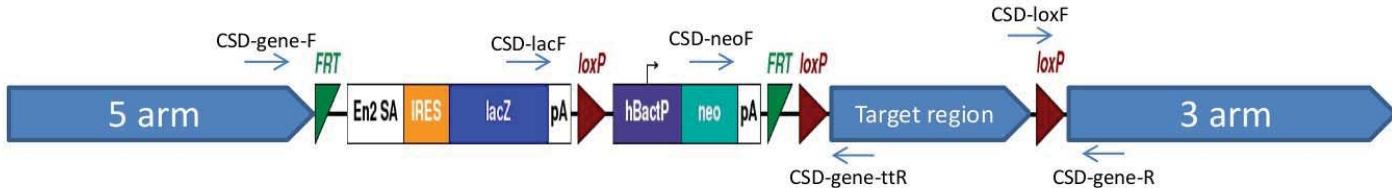

Cassette Primers

CSD-lacF: GCTACCATTACCAAGTTGGTCTGGTGTC
CSD-neoF: GGGATCTCATGCTGGAGTTCTTCG
CSD-loxF: GAGATGGCGAACGCAATTAATG

Gene Specific Primers

CSD-R: GAAGCAAAGAAACAGGAAAACACTGG
CSD-ttR: CCTGACCTTAACCTAGGGGACTCCG
CSD-F: ACGAAATTAAAAACTTAAGCCTGTCC

| Genotype | Forward Primer | Reverse Primer | Amplicon size (bp) |
|---------------|----------------|----------------|--------------------|
| Floxed | CSD-loxF | CSD-R | 264 |
| PreCre | CSD-neoF | CSD-ttR | 584 |
| PostCre | CSD-lacF | CSD-R | 562 |
| Wildtype | CSD-F | CSD-ttR | 326 |
| PostFlp | CSD-F | CSD-ttR | 513 |
| PostFlp & Cre | CSD-F | CSD--R | 540 |

Please note, these primers are auto-designed and may not have been verified by the repository, and as such may require optimization or redesign by your facility.

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the Floxed or PreCre primers may be used to identify the mutant mice. The Floxed primers test for the distal LoxP site. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line (without any FLP recombination). The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The PostFlp & Cre primers should be used if mutant mice were crossed with a Flp recombinase line and then a Cre recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).