

MMRRC Stock #: 047074-UCD

C57BL/6N-Atm1Brd Ckap69tm1a(KOMP)Wtsi/  
Mmucd

**Mouse PCR protocol**

Design ID: 44774

selection cassette: L1L2 Bact P


Suggested DNA Prep: DNeasy® Tissue Kit

Reagents	1X (uL)
Water (biology grade)	10.725
betain 5M (Sigma)	6.5
dmsO (Sigma)	0.325
10X Buffer w/o MgCl <sub>2</sub> (AB)	2.5
25 mM MgCl <sub>2</sub> (AB)	1.75
10 mM dNTPs (invitrogen)	0.5
Primers (20uM each)	0.5
taq 5U/uL (amplitaq, AB)	0.2
total cocktail	23
template	2
reaction volume	25

**Cycling Parameters**

Temp C	Time	
94	5 min	
94	15 sec	
65	30 sec	10X (decrease 1C/cycle)
72	40 sec	
94	15 sec	
55	30 sec	30X
72	40 sec	
72	5 min	
4	finished	

**Primer strategy**

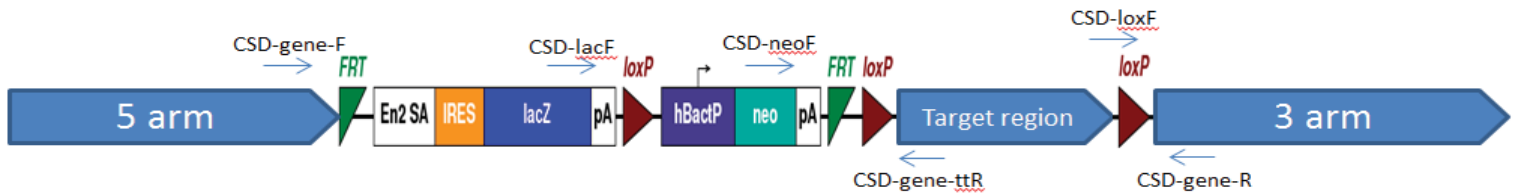

**Cassette Primers**

CSD-lacF: gctaccattaccagttggtctggtgctc

CSD-neoF: gggatctcatgctggagttctctcg

CSD-loxR: gagatggcgcaacgcaattaatg

**Gene Specific Primers**

CSD-A330021E22Rik-R: TTTCCTTGGGTGCTCTCTCACTCCC

CSD-A330021E22Rik-ttR: CGCTGCAACTGGAATCAGAAGTAGG

CSD-A330021E22Rik-F2: caaacagcacaggagattcaaatgg

Genotypes	Primer pairs	amplicon size (bp)
Floxed	loxR & -R	398
PreCre	neoF & -ttR	581
PostCre	lacF & -R	696
Wildtype	-F2 & -ttR	364
PostFlp	-F & -ttR	
PostFlp & Cre	-F & -R	

Please note, these primers are auto-designed and may not have been verified by the repository, and as such may require optimization or redesign by your facility.

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the Floxed or PreCre primers may be used to identify the mutant mice. The Floxed primers test for the distal LoxP site. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line (without any FLP recombination). The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The PostFlp & Cre primers should be used if mutant mice were crossed with a Flp recombinase line and then a Cre recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).