

KOMP PCR Design

Mouse PCR Protocol (version 1 PK)

Design ID: 208473
Project ID: CSD72177
Selection Cassette: L1L2_Bact_P

Suggested DNA Prep: DNeasy®Tissue Kit

Reagent	1X (μL)
water (biological grade)	10.725
betain 5M (Sigma)	6.5
DMSO (Sigma)	0.325
10X buffer w/o MgCl ₂ (AB)	2.5
25 mM MgCl ₂ (AB)	1.75
10 mM dNTPs (Invitrogen)	0.5
primers (20 μM each)	0.5
Taq 5U/μL (AmpliTaq, AB)	0.2
total cocktail	23
template	2
reaction volume	25




MMRRC Stock #: 048162-UCD

C57BL/6N-Atm1Brd Igflr1tm1a(KOMP)Mbp/
TcpMmucd

Cycling Parameters	
Temperature °C	Time
94	5 min
94	15 sec
65	30 sec
72	40 sec
10X (decrease 1°C/cycle)	
94	15 sec
55	30 sec
72	40 sec
30X	
72	5 min
4	finished

Primer Strategy

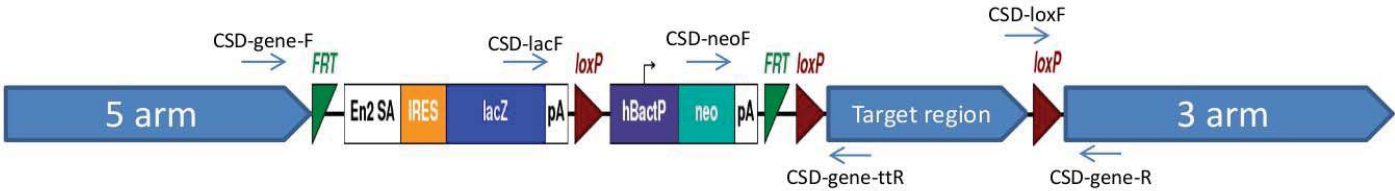

Cassette Primers

CSD-lacF: GCTACCATTACCAGTTGGTCTGGTGTC
CSD-neoF: GGGATCTCATGCTGGAGTTCTTCG
CSD-loxF: GAGATGGCGCAACGCAATTAATG

Gene Specific Primers

CSD-R: GCCAGGACTGTTGACATATTTGCATATC
CSD-ttR: TAAAGCCTACCTAAGTCCATTCTCCC
CSD-F: actggaactctgacaagaggtgctg

Genotype	Forward Primer	Reverse Primer	Amplicon size (bp)
Floxed	CSD-loxF	CSD-R	341
PreCre	CSD-neoF	CSD-ttR	598
PostCre	CSD-lacF	CSD-R	639
Wildtype	CSD-F	CSD-ttR	311
PostFlp	CSD-F	CSD-ttR	529
PostFlp & Cre	CSD-F	CSD-R	619

Please note, these primers are auto-designed and may not have been verified by the repository, and as such may require optimization or redesign by your facility.

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the Floxed or PreCre primers may be used to identify the mutant mice. The Floxed primers test for the distal LoxP site. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line (without any FLP recombination). The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The PostFlp & Cre primers should be used if mutant mice were crossed with a Flp recombinase line and then a Cre recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).