

KOMP PCR Design

Mouse PCR Protocol (version 1JE)

Design ID: 371847
Project ID: CSD93064
Selection Cassette L1L2_Bact_P



MMRRC Stock #: 048941-UCD

C57BL/6N-Atm1Brd P3h1tm1a(EUCOMM)Wtsi/
MbpMmucd


Suggested DNA Prep: DNeasy®Tissue Kit

Reagent	1X (μL)
water (biological grade)	10.725
betain 5M (Sigma)	6.5
DMSO (Sigma)	0.325
10X buffer w/o MgCl ₂ (AB)	2.5
25 mM MgCl ₂ (AB)	1.75
10 mM dNTPs (Invitrogen)	0.5
primers (20 μM each)	0.5
Taq 5U/μL (AmpliTaq, AB)	0.2
total cocktail	23
template	2
reaction volume	25

Cycling Parameters

Temperature °C	Time	
94	5 min	
94	15 sec	10X (decrease 1°C/cycle)
65	30 sec	
72	40 sec	
94	15 sec	30X
55	30 sec	
72	40 sec	
72	5 min	
4	finished	

Primer Strategy

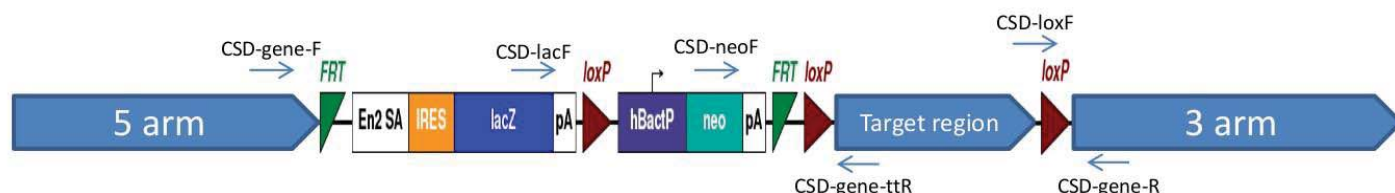

Cassette Primers

CSD-lacF: GCTACCATTACCAAGTTGGTCTGGTGTC
CSD-NeoF: GGGATCTCATGCTGGAGTTCTTCG
CSD-lacF: GAGATGGCGCAACGCAATTAATG

Gene Specific Primers

CSD-R: TGTGGGAGGAACTGTGGCTTG
CSD-ttR: AGTAAACAAGCTGATTTGCTGAGG
CSD-F: GCATGTCTAGCAGCTATGTCTGTGG

Genotype	Forward Primer	Reverse Primer	Amplicon size (bp)
Floxed	CSD-lacF	CSD-R	258
PreCre	CSD-NeoF	CSD-ttR	618
PostCre	CSD-lacF	CSD-R	556
Wildtype	CSD-F	CSD-ttR	311
PostFlp	CSD-F	CSD-ttR	506
PostFlp & Cre	CSD-F	CSD-R	493

Please note, these primers are auto-designed and may not have been verified by the repository, and as such may require optimization or redesign by your facility.

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the Floxed or PreCre primers may be used to identify the mutant mice. The Floxed primers test for the distal LoxP site. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line (without any FLP recombination). The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The PostFlp & Cre primers should be used if mutant mice were crossed with a Flp recombinase line and then a Cre recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).