

# KOMP PCR Design

Mouse PCR Protocol (version 1)

Design ID: 278732

Project ID: CSD102779

Selection Cassette: L1L2\_Bact\_P



MMRRC Stock #: 049192-UCD

C57BL/6N-Atm1Brd Preptm1a(KOMP)Wtsi/  
JMmucd


## Suggested DNA Prep: DNeasy® Tissue Kit

Reagent	1X (μL)
water (biological grade)	10.725
betain 5M (Sigma)	6.5
DMSO (Sigma)	0.325
10X buffer w/o MgCl <sub>2</sub> (AB)	2.5
25 mM MgCl <sub>2</sub> (AB)	1.75
10 mM dNTPs (Invitrogen)	0.5
primers (20 μM each)	0.5
Taq 5U/μL (AmpliTaq, AB)	0.2
total cocktail	23
template	2
reaction volume	25

## Cycling Parameters

Temperature °C	Time	
94	5 min	
94	15 sec	10X (decrease 1°C/cycle)
65	30 sec	
72	40 sec	
94	15 sec	30X
55	30 sec	
72	40 sec	
72	5 min	
4	finished	

## Primer Strategy

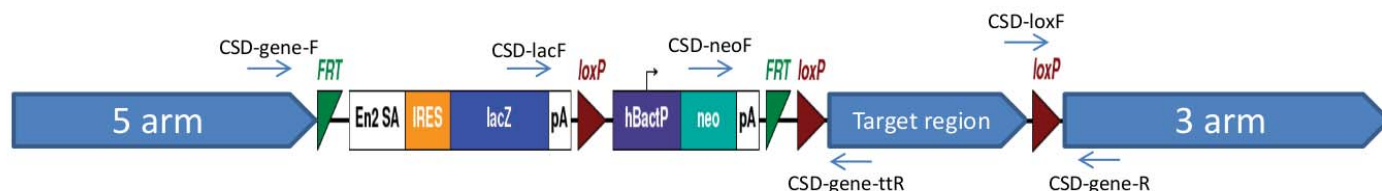

## Cassette Primers

CSD-lacF: GCTACCATTACCAGTTGGTCTGGTGTC  
 CSD-neoF: GGGATCTCATGCTGGAGTTCTTCG  
 CSD-loxF: GAGATGGCGCAACGCAATTAATG

## Gene Specific Primers

CSD-Prep-R: AACACGGCCAACCTATTTCAAGAGC  
 CSD-Prep-ttR: CTGTCAGTGATGGAGTGTGACCTGG  
 CSD-Prep-F: TGTCATGCATTCTGTGTGCTCTTGG

## Genotype Forward Primer Reverse Primer Amplicon size (bp)

Genotype	Forward Primer	Reverse Primer	Amplicon size (bp)
Floxed	CSD-loxF	CSD-Prep-R	325
PreCre	CSD-neoF	CSD-Prep-ttR	646
PostCre	CSD-lacF	CSD-Prep-R	623
Wildtype	CSD-Prep-F	CSD-Prep-ttR	580
PostFlp	CSD-Prep-F	CSD-Prep-ttR	775
PostFlp & Cre	CSD-Prep-F	CSD-Prep-R	747

Please note, these primers are auto-designed and may not have been verified by the repository, and as such may require optimization or redesign by your facility.

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the Floxed or PreCre primers may be used to identify the mutant mice. The Floxed primers test for the distal LoxP site. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line (without any FLP recombination). The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The PostFlp & Cre primers should be used if mutant mice were crossed with a Flp recombinase line and then a Cre recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).