


Mouse PCR Protocol (version 1)

Selection Cassette: L1L2_Bact_P

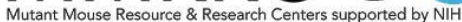

C57BL/6N-Rbmxtm2b(KOMP)Wtsi/Mmucd

Reagent	1X (μL)
---------	----------

water (biological grade)	10.725
betain 5M (Sigma)	6.5
DMSO (Sigma)	0.325
10X buffer w/o MgCl ₂ (AB)	2.5
25 mM MgCl ₂ (AB)	1.75
10 mM dNTPs (Invitrogen)	0.5
primers (20 μM each)	0.5
Taq 5U/μL (AmpliTaq, AB)	0.2
total cocktail	23
template	2
reaction volume	25

Temperature °C	Time
----------------	------

94	5 min	
94	15 sec	
65	30 sec	10X (decrease 1°C/cycle)
72	40 sec	
94	15 sec	
55	30 sec	30X
72	40 sec	
72	5 min	
4	finished	

CSD-lacF: GCTACCATTACCAGTTGGTCTGGTGTCT
CSD-neoF: GGGATCTCATGCTGGAGTTCTTCG
CSD-loxF: GAGATGGCGCAACGCAATTAATG

CSD-RbmX-R:	GTCTCTCGTTTTACTGGCAGAGGTCC
CSD-RbmX-ttR:	CTCTGTATAGCTGTGGCTGTCCTGG
CSD-RbmX-F:	CGGTGTGTCCAGAAATAGACCCATGC

Genotype	Forward Primer	Reverse Primer	Amplicon size (bp)
Floxed	CSD-loxF	CSD-RbmX-R	397
PreCre	CSD-neoF	CSD-RbmX-ttR	530
PostCre	CSD-lacF	CSD-RbmX-R	695
Wildtype	CSD-RbmX-F	CSD-RbmX-ttR	426
PostFlp	CSD-RbmX-F	CSD-RbmX-ttR	570
PostFlp & Cre	CSD-RbmX-F	CSD-RbmX-R	730

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the Floxed or PreCre primers may be used to identify the mutant mice. The Floxed primers test for the distal LoxP site. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line (without any FLP recombination). The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The PostFlp & Cre primers should be used if mutant mice were crossed with a Flp recombinase line and then a Cre recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).