

Mouse PCR Protocol (version 1BW)

Design ID: 210835

Project ID: CSD69150

Selection Cassette: L1L2_Bad_P_



MMRRC Stock #: 049425-UCD

C57BL/6N-Atm1Brd Rgcctm1(KOMP)Wtsi/


MbpMmucd

Suggested DNA Prep: DNeasy®Tissue Kit

Reagent	1X (μ L)
water (biological grade)	10.725
betain 5M (Sigma)	6.5
DMSO (Sigma)	0.325
10X buffer w/o MgCl ₂ (AB)	2.5
25 mM MgCl ₂ (AB)	1.75
10 mM dNTPs (Invitrogen)	0.5
primers (20 μ M each)	0.5
Taq 5U/ μ L (AmpliTaq, AB)	0.2
total cocktail	23
template	2
reaction volume	25

Cycling Parameters		
Temperature °C	Time	
94	5 min	
94	15 sec	
65	30 sec	10X (decrease 1°C/cycle)
72	40 sec	
94	15 sec	
55	30 sec	30X
72	40 sec	
72	5 min	
4	finished	

Primer Strategy

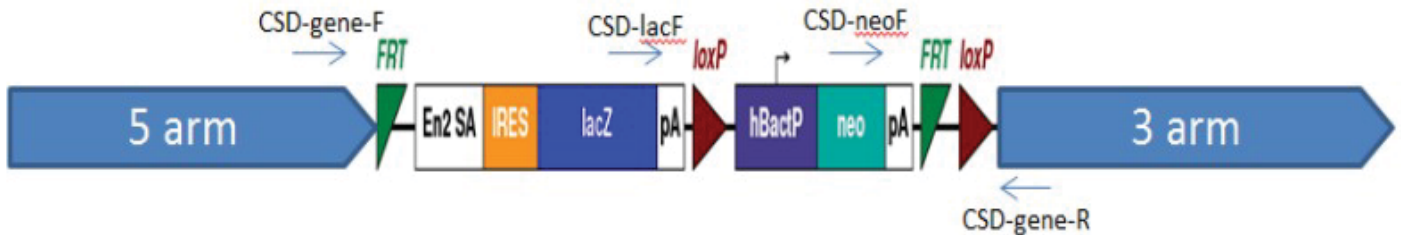

Cassette Primers

CSD-lacF: GCTACCATTACCAGTTGGTCTGGTGTC

CSD-neoF: GGGATCTCATGCTGGAGTTCTTCG

Gene Specific Primers

CSD-R: GGGAAAGCCTAGTTGAGGATACGG

CSD-F: TCTGGTGACAGGCGCTTTAGAGG

CSD-wtR: AACTCGCACAGCACGTCGGTC

Genotype	Forward Primer	Reverse Primer	Amplicon size (bp)
PreCre	CSD-neoF	CSD-R	525
PostCre	CSD-lacF	CSD-R	476
Wildtype	CSD-F	CSD-wtR	194
PostFlp	CSD-F	CSD-R	346

Please note, these primers are auto-designed and may not have been verified by the repository, and as such may require optimization or redesign by your facility.

We recommend running primers singleplex. For screening of pups prior to any Flp or Cre recombination, the PreCre primers may be used to identify the mutant mice. The PostCre primers should be used if mutant mice were crossed with a Cre recombinase line. The PostFlp primers should be used if mutant mice were crossed with a Flp recombinase line. The wildtype primers should be used for zygosity testing of all mutant mice (pre or post recombination).