

GENOTYPING BY PCR PROTOCOL

MUTANT MOUSE RESOURCE & RESEARCH CENTER: UC DAVIS

mmrrc@ucdavis.edu

530-754-MMRRC

Protocol Name: B6.Cg-Tbr1tm2.1Jlr/Mmucd **MMRRC: 065429-UCD**

Protocol:

Reagent/Constituent	Volume (μL)
Water	5.6
GoTaq® G2 Colorless Master Mix, 2X	7.5
Primer 1. (stock concentration is 20μM)	0.45
Primer 2. (stock concentration is 20μM)	0.45
DNA (example) extracted w/ "Qiagen DNeasy columns or other similar silica based kits"	1.0
TOTAL VOLUME	
15	

Comments on protocol:

- Protocol may work with other DNA extraction methods.
- Use Touch-Down cycling protocol-first 10 cycles anneal at 65°C decreasing in temperature by 1.0°C; next 30 cycles anneal at 55°C.
- The mutant PCR is a general LacZ PCR. The wild type is specific for this strain.

Strategy:


Steps	Temp (°C)	Time (m:ss)	# of Cycles
1. Initiation/Melting HOT START? ☐	94	2:00	1x
2. Denaturation	94	0:10	
3. Annealing steps 2-3-4 cycle in sequence	65 (↓1°C/cycle)	0:30	10x
4. Elongation	68	2:00	
5. Denaturation	94	0:15	
6. Annealing steps 5-6-7 cycle in sequence	55	0:30	25x
7. Elongation	68	2:00 (↑20sec/cycle)	

Primers:

Electrophoresis Protocol:

Name	Nucleotide Sequence (5' - 3')	Argarose: 1.5% V: 90		
1. 65429-Tbr1-loxF	GACACACACCCTTCTTCAGTTTACAGC	Estimated Running Time: 90 min.		
2. 65429-Tbr1-loxR	CAAGCCCGACTGCCAATGTTCTG	Primer Combination	Band (bp)	Genotype
		1 & 2	400	mutant
		1 & 2	280	wildtype

GENOTYPING BY PCR PROTOCOL MUTANT MOUSE RESOURCE & RESEARCH CENTER: UC DAVIS mmrrc@ucdavis.edu 530-754-MMRRC

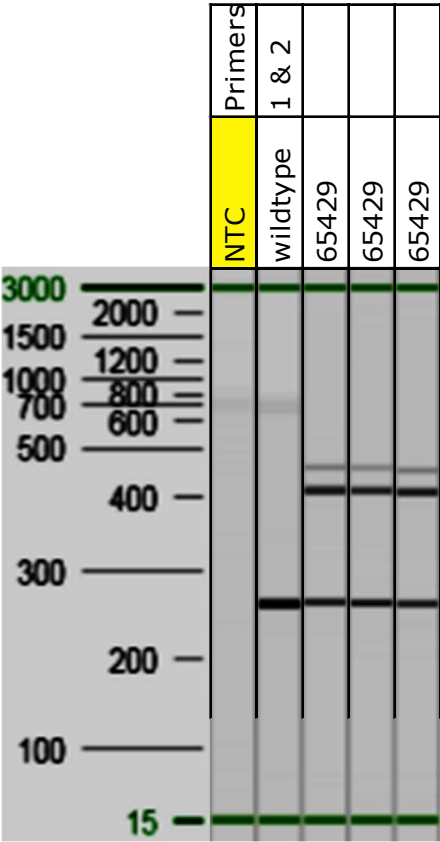

Sample ID	Assay	Well
ntc		A1
wt		B1
10466		A2
2		B2
3		A3
4		B3
5		A4
6		B4
7		A5
8		B5
9		A6
ntc		B6
wt		A7
32358		B7
3		A8
4		B8
5		A9
6		B9
7		A10
8		B10
9		A11
10		B11
11		A12
12		B12
13		C1
14		D1
15		C2
ntc		D2
wt		C3
31665		D3
3		C4
4		D4
5		C5
6		D5
7		C6
8		D6

PCR protocol