

GENOTYPING PROTOCOL

MUTANT MOUSE RESOURCE & RESEARCH CENTER: UC DAVIS

mmrrc@ucdavis.edu

530-754-MMRRC

Protocol Name: CR10520 Vcan EXDEL

Protocol: GoTaq® G2 Colorless Master Mix(Promega)

Reagent/Constituent	Volume (μL)
Water	4.5
GoTaq® G2 Colorless Master Mix,2X	7.5
Primer 1. (stock concentration is 20μM) comF	0.5
Primer 2. (stock concentration is 20μM) wtR	0.5
Primer 3. (stock concentration is 20μM) mutR	0.5
DNA (example) extracted w/ "Qiagen DNeasy columns or other similar silica based kits"	1.5
TOTAL VOLUME OF REACTION:	15.00 μL

Comments on protocol:

- Protocol may work with other DNA extraction methods.

Strategy:

Steps	Temp (°C)	Time (m:ss)	# of Cycles
1. Initiation/Melting HOT START? ☐	94	2:00	1x
2. Denaturation	94	0:10	
3. Annealing steps 2-3-4 cycle in sequence	65 (↓1°C/cycle)	0:30	10x
4. Elongation	68	2:00	
5. Denaturation	94	0:15	
6. Annealing steps 5-6-7 cycle in sequence	55	0:30	25x
7. Elongation	68	2:00 (↑20sec/cycle)	
8. Finish	4	∞	n/a


Primers:

Name	Nucleotide Sequence (5' - 3')	Agarose: 1.5% V: 90
1. CR_Vcan-5IVF (comF)	TTCTACAATAGCATTCCTCTGGGCA	Estimated Running 90 min.
2. CR_Vcan-5IVR (wtR)	CCGCTCTCATCTTCCTCACTGTGAT	Primer Combination Band (bp) Genotype
3. CR_Vcan-3IVR (mutR)	CATGGTATGTCTGCACATCCTTGG	1 & 2, 1 & 3 1225, 7019 wildtype
		1 & 3 1161 mutant


Electrophoresis Protocol:

Allele Description: Exon 8 [ENSMUSE00000680689](#) and flanking splicing regions were constitutively deleted from the Vcan gene [ENSMUSG0000021614.16](#) using CRISPR Cas9 gene editing technology in mouse zygotes. Subsequent founders were backcrossed to C57BL6/N to produce sequence confirmed heterozygous animals.

WT:

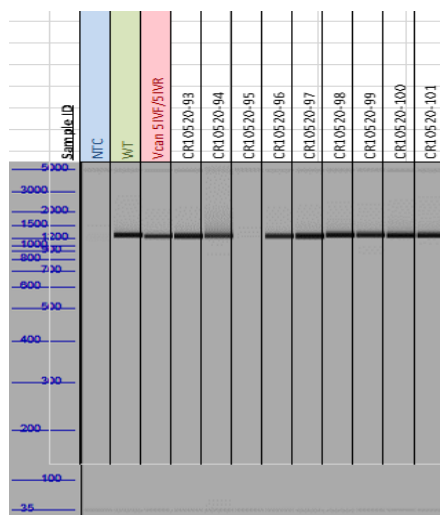

Mut:

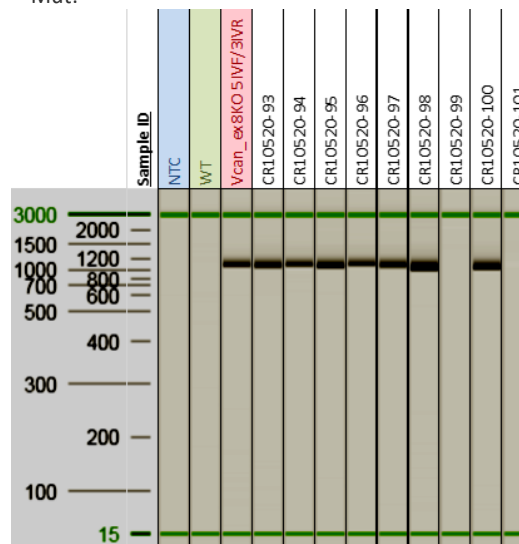

WT = 99 and 101; Het = 92-94, 96-98, 100; Hom = 95