

GENOTYPING PROTOCOL

MUTANT MOUSE RESOURCE & RESEARCH CENTER: UC DAVIS

mmrrc@ucdavis.edu

530-754-MMRRC

Protocol Name: CR10769 Id3 exdel

Protocol: GoTaq® G2 Colorless Master Mix(Promega)

Reagent/Constituent	Volume (μL)
Water	4.5
GoTaq® G2 Colorless Master Mix,2X	7.5
Primer 1. (stock concentration is 20μM) comF	0.5
Primer 2. (stock concentration is 20μM) wtR	0.5
Primer 3. (stock concentration is 20μM) mutR	0.5
DNA (example) extracted w/ "Qiagen DNeasy columns or other similar silica based kits"	1.5
TOTAL VOLUME OF REACTION:	15.00 μL

Comments on protocol:

- Protocol may work with other DNA extraction methods.

Strategy:

Steps	Temp (°C)	Time (m:ss)	# of Cycles
1. Initiation/Melting HOT START? ☐	94	2:00	1x
2. Denaturation	94	0:10	
3. Annealing steps 2-3-4 cycle in sequence	65 (↓1°C/cycle)	0:30	10x
4. Elongation	68	2:00	
5. Denaturation	94	0:15	
6. Annealing steps 5-6-7 cycle in sequence	55	0:30	25x
7. Elongation	68	2:00 (↑20sec/cycle)	
8. Finish	4	∞	n/a


Primers:

Electrophoresis Protocol:

Name	Nucleotide Sequence (5' - 3')	Agarose: 1.5%	V: 90
1. CR Id3 comF	GTCTGGAGGTCAGACGAGCAGCA	Estimated Running Time: 90 min.	
2. CR Id3 wtR	TCGAGGATGTAGTCTATGACACGCTG	Primer Combination	Band (bp)
3. CR Id3 mutR	CCAGAAGAACAGCTCTTATGCTGCC	1 & 2, 1 & 3	789, 1902
		1 & 3	755
			wildtype
			mutant

Allele Description: The coding region of Exon 1 ENSMUSE00000795245, the entirety of Exon 2 ENSMUSE00000181163, and flanking splicing regions were constitutively deleted from the Id3 gene [ENSMUST00000008016.2](#) using CRISPR Cas9 gene editing technology in mouse zygotes. Subsequent founders were backcrossed to C57BL6/N to produce sequence confirmed heterozygous animals.

*wtR primer untested (ePCR verified)

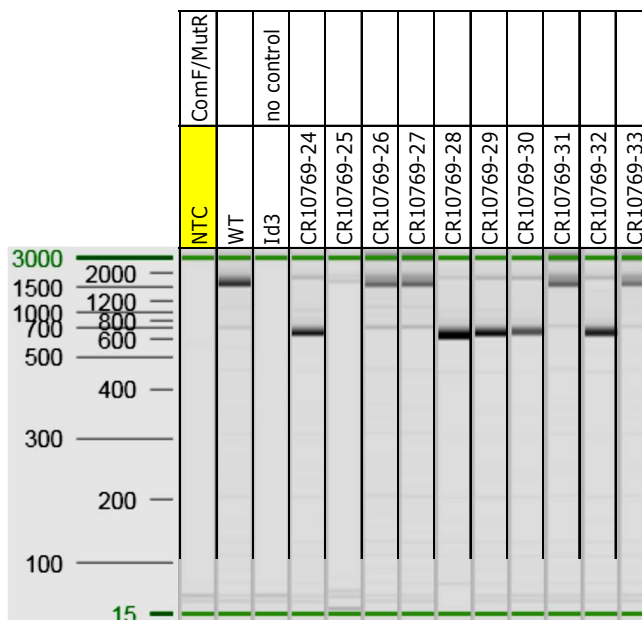

PCR protocol developed by MMRRC at University of California, Davis