

# GENOTYPING BY PCR PROTOCOL

## MUTANT MOUSE RESOURCE & RESEARCH CENTER: UC DAVIS

[mmrrc@ucdavis.edu](mailto:mmrrc@ucdavis.edu)  
530-754-MMRRC

**STRAIN NAME:** B6.Cg-Tg(CD2-Tcra,Cd4-Tcrb)4CKang/Mmucd **MMRRC:** 037367-UCD

**Protocol:**

Reagent/Constituent	Volume (µL)
Water	11.275
10x Buffer	2.5
MgCl <sub>2</sub> (stock concentration is 25mM)	1.7
Betaine (stock concentration is 5M) <i>Optional</i>	6.5
dNTPs (stock concentration is 10mM)	0.5
DMSO <i>Optional</i>	0.325
Primer 1. (stock concentration is 20µM)	0.5
Primer 2. (stock concentration is 20µM)	0.5
Taq Polymerase 5Units/µL	0.2
DNA (example) extracted w/ "Qiagen DNeasy columns or other similar silica based kits"	1.0
<b>TOTAL VOLUME OF REACTION:</b>	<b>25.000 µL</b>

**Comments on protocol:**

- Protocol may work with other DNA extraction methods.
- Use Touch-Down cycling protocol-first 10 cycles anneal at 65°C decreasing in temperature by 1.0°C; next 30 cycles anneal at 55°C.
- Betaine and DMSO have been standardized due to high GC content. Protocol may be tested without. Also, may adjust


**Strategy:**

Steps	Temp (°C)	Time (m:ss)	# of Cycles
1. Initiation/Melting <span style="float: right;">HOT START? <input type="checkbox"/></span>	94	5:00	1
2. Denaturation	94	0:15	
3. Annealing <span style="float: right;">steps 2-3-4 cycle in sequence</span>	65 to 55 (↓1°C/cycle)	0:30	<b>40x</b>
4. Elongation	72	0:40	
5. Amplification	72	5:00	1
6. Finish	15	∞	n/a

**Primers:**

**Electrophoresis Protocol**

Name	Nucleotide Sequence (5' - 3')	Agarose: 1.5%	V: 90
1. 37367 tcrA-tgF	CTGACACCTGCTCAGTTCCTTGTC	Estimated Running Time: 90 min.	
2. 37367 tcrA-tgR	GGTGCAGTTCAGCATCACAGG	<b>Primer Combination</b>	<b>Band (bp)</b>
3. 37367_tcrB-tgF	AGAGGATGGTGGCAGACAAGACC	1 and 2	122
4. 37367_tcrB-tgR	GGCTTCAGAGGAGGACAAGTGG	3 and 4	141
			<b>Genotype</b>
			Tcra Tg
			Tcrb Tg

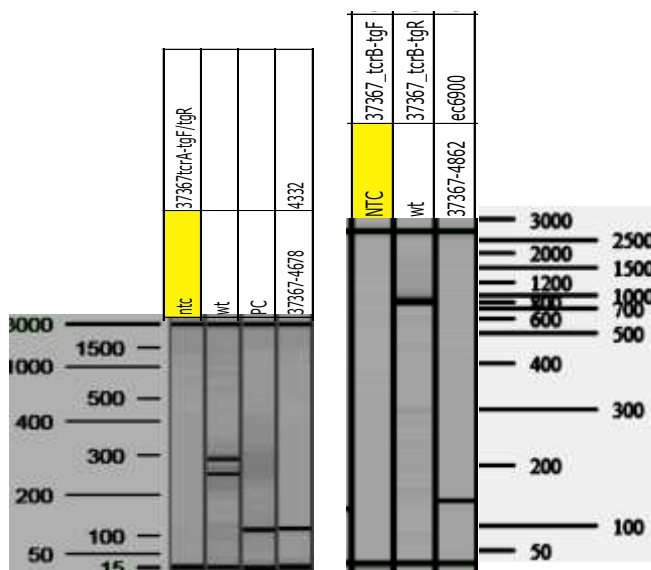

**Please note:** *Tcra* Primers will pick up 300/260 bp WT background and *Tcrb* will pick 750bp